

FAPAS for allergens in Cake Mix

Results from *GlutenTox PRO* AOAC-RI Certified (Biomedal)

Allergen	Round	Identification	Date	Lab identification number	Type of matrix	Biomedal's results	Consensus assigned value	Participants (number)	Max-Min values
Gluten	fapas	27204	Aug 2017	134	Cake Mix A (TMA)	<u>>20 ppm</u> <u><40 ppm</u>	<u>24 ppm</u>	143	49.6-11.8 ppm
Gluten	fapas	27204	Aug 2017	134	Cake Mix B (TMB)	<u>>40 ppm</u>	<u>41.3 ppm</u>	143	107-17.3 ppm

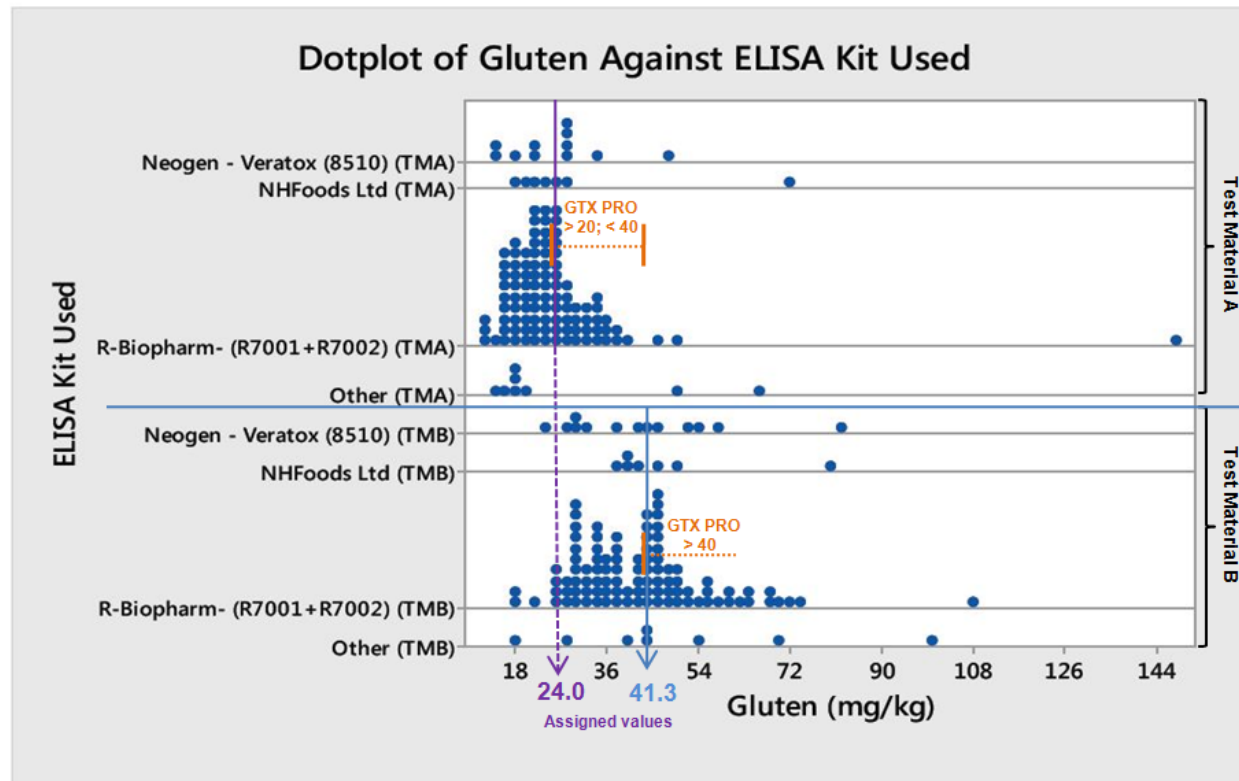


Figure 25: Dot Plot of Gluten Data in Cake Mix by Kit Type

Figure 25. Semiquantitative Results of *GlutenTox PRO* with test material A (TMA) and B (TMB)